


The Set-Up Procedure



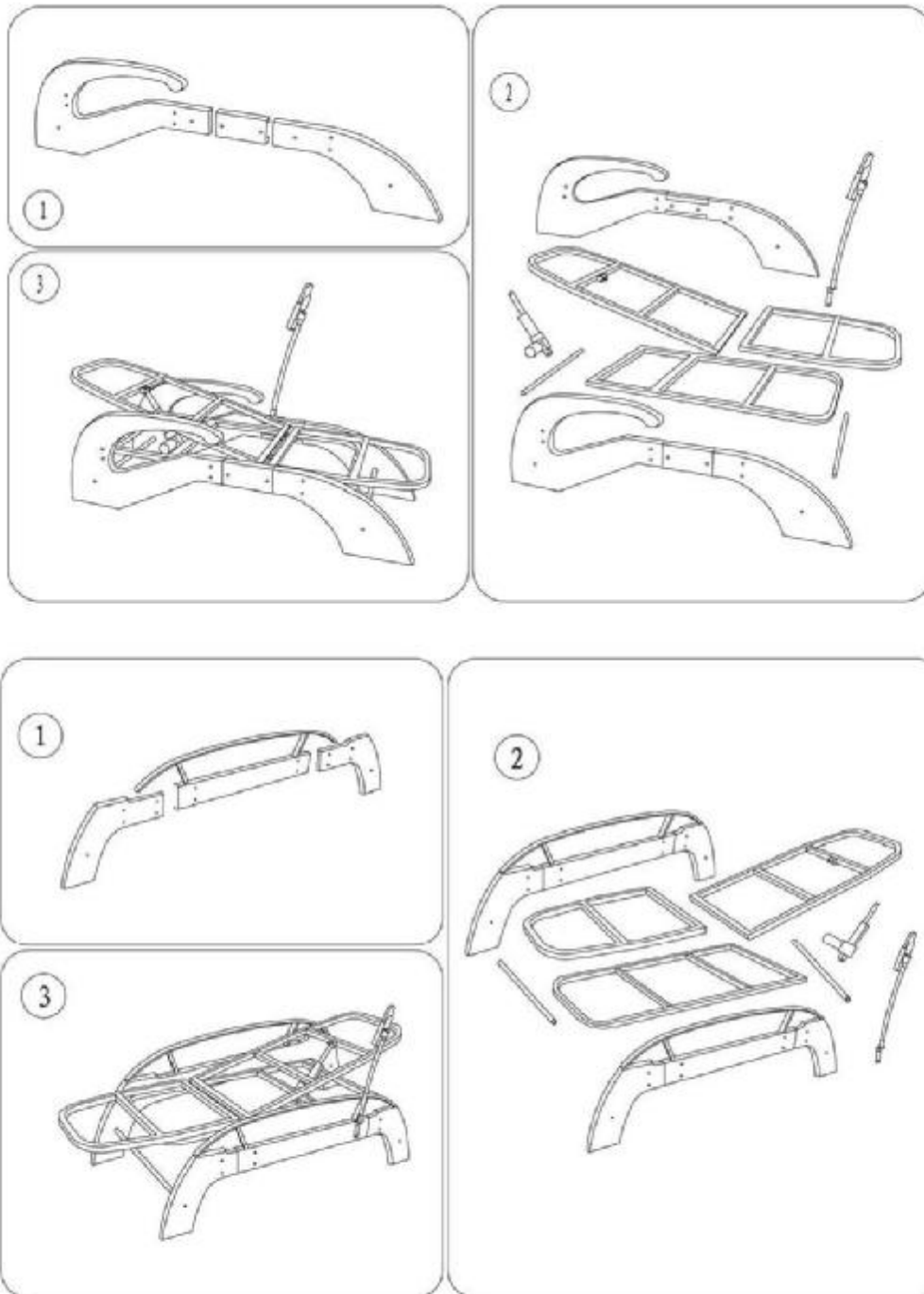

- The title page and the back cover are the identical function product, in the use operation also same, only cosmetics difference.

Specifications are subject to change without notice or obligation on the part of manufacture

User Manual

Thank you for using our product,
Please read this manual carefully before using.



CE
RoHS

Jade Thermal Massage Bed

Table Of Contents

Product Introduction	02
General Product Information	
Main Features	
Messages For Users	
Indications	
Contraindications	
Precautions	
Inspection Before Using The Product	
Inspection After Using The Product	
Product Exploded View	
Product Drawing	
Remote Control	
How To Use Our Product	07
Using The Product	
Instructions For Using LED Controller	
Tension Adjustment	
Installation Of The Product In Your Home	11
Points To Remember Before Installation	
Installation Instruction	
Maintenance & Cautions	11
Safety Practice	
Safety	
Care And Storage	
Remote Control Operation & Error Messages Explanations	
Troubleshooting & Service	
Self Troubleshooting	
Service	
Specifications	
Additional information	14
The Set-Up Procedure	15

Product Introduction

General Product Information

The manufacturing of this product (Proleader PLD-6018 series Jade Thermal Massage Bed) has been approved by Chinese SFDA (State Food and Drug Administration). It is designed to target your body's vital pressure points along the spine with a precise, relaxing massage to provide you with the combined therapeutic benefits of massage, chiropractic, acupressure and Far-Infrared Heat. It relieves stress and tension in the muscles, improve blood circulation, and increase energy.

Main Features

1. Four or six groups of premium jade internal projector for maximum effect
2. Two external jade projectors with 5/9 balls or 5/13 balls to be applied on the parts of body unreachable by the bed
3. Back massage mattress can be adjusted to 45 degree angle to suit for your personal needs
4. Heated jade roller reach temperatures from 40C to 70C to effectively transfer heat to your needed discomfort spots on your back
5. Neck jade roller with linear rail to provide extra treatment for the head and neck discomfort
6. Tourmaline stone has five major benefit for human body:
 - a. Release of 4-14um high energy far infrared ray;
 - b. Release of negative ion permanently;
 - c. Generate 0.06ma bio-current;
 - d. Rich in minerals and elements?
 - e. Erase bad smells and kills bad bacteria
7. Chiropractic, acupressure and infrared therapy

Messages for Users

Indications

Used as supplementary instrument treatment for lumbar spondylosis (chronically strain of lumbar muscles, slipped lumbar disc, back pain) and cervical spine problems.

Contraindications

People with follow illness are strictly prohibited to use this product:

1. In pregnancy
2. With tumor
3. Severe high blood pressure
4. Acute suppurative inflammation
5. Hemorrhagic tendency
6. Artificial cardiac pacemaker
7. Bone fracture

Precautions

1. Cotton clothes are recommended during therapy. Please remove the belt and tie before therapy.
2. Please do not place external projectors directly on the skin to avoid the possible burn in case of excessive heat resulted by improper setting of the temperature.
3. Please do not place external projectors directly over your face to avoid the contact with eyes.
4. If you are sensitive to heat, please lower the temperature to proper level before using the product.
5. Please follow the instructions and do not use other non-Proleader products while using Proleader products.
6. Please do not go to sleep during the program.
7. This device can not replace any other medical treatments/surgeries suggested or prescribed by your doctor for curing/healing purpose. If you symptom persists or worsen, please seek for medical help.

Inspection Before Using The Product

1. Make sure your power cable is plugged in a working power outlet;
2. Check if external projectors are plugged in properly with the bed;


3. Check if main power switch and remote control work properly;
4. Check if temperature display is correct on the remote control;
5. Check if all the lamps are working properly;
6. Please make sure the projector's temperature is around 45C–55C;
7. Check if any abnormality about projectors







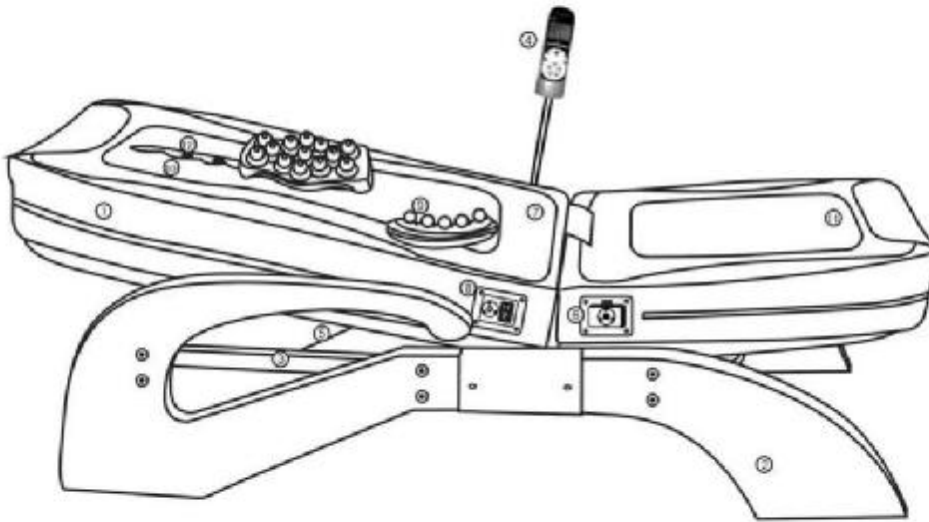


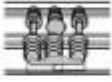
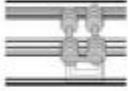

Inspection After Using The Product

1. Check if temperature display is turned off;
2. Check if all lamps are shut off;
3. Check if power switch is at off position;
4. Check if power cable is unplugged from power outlet;

Product Exploded View

Product Drawing



 1. Mattress	 2. Supporting leg	 3. Bed frame	 4. Remote control
 5. Internal projector	 6. Leg outlet	 7. Panel for remote control and external projectors	 8. Power panel
 9. External handheld projectors	 10. Internal projector	 11. Leg part	 12. Magnetic patch


Remote Control



How to use our product

Using The Product

1. See the picture below. Lay with your belly up on the massage bed. Keep your head a fist down the top of the massage bed. Rest your legs and arms on the mattress.
2. When using Automatic modes, please place the external projectors on your desired spots. (Refer to below description for the safety practice)
3. You may continuously use Automatic or Manual mode to reach the desired results.

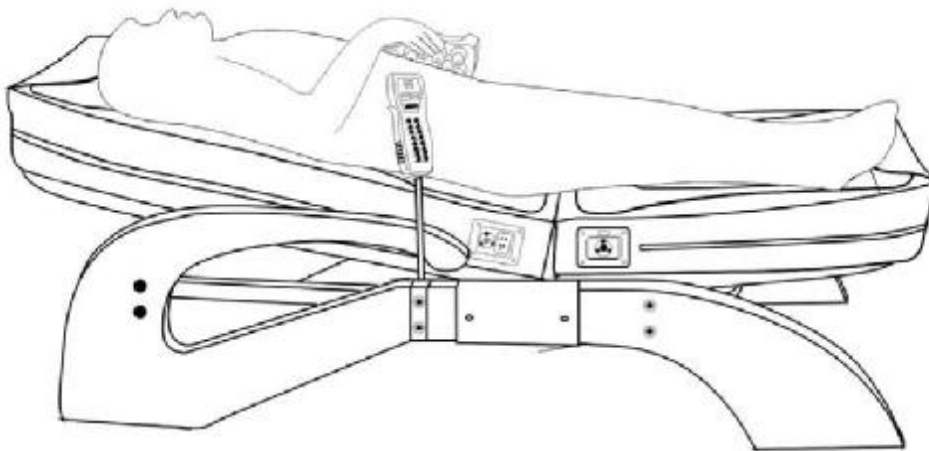

Instructions For Using LED Controller

Press Power/Preheating/Automatic Key

Back projector, handheld projector and infrared start a 10 minutes of preheating process.LED window displays the current temperatures.After preheating, press Power/Preheating/Automatic Key again, massage bed enters N1 default program Automatic mode.

Press Back Automatic Key to choose among modes of N1, N2, N3 and N0. N1, N2, N3 are the automatic modes and N0 is manual mode.


During N0 manual mode, use the left and right arrow Back Manual Key to move the back projectors to desired spots.



Press Leg Automatic Key to choose between automatic and manual modes. Use left and right arrow Leg Manual Key to move the leg projectors to desired spots in manual mode.

Use left and right Back Angle Adjustment Key to rise the back stand to desired angle.

Temperature setting:

Use the temperature dial to set the temperatures for back projector, infrared and handheld projector. To set temperature for the back projector, press Back Projector Key once. On the LED screen, the back projector temperature will blink.Press Temperature Adjustment key to increase or decrease of the temperature. Press + or - once will cause temperature to rise or drop in 5°C. (Temperature range is 40°C - 70°C) Use the same procedure to set temperature for Infrared and handheld projector.

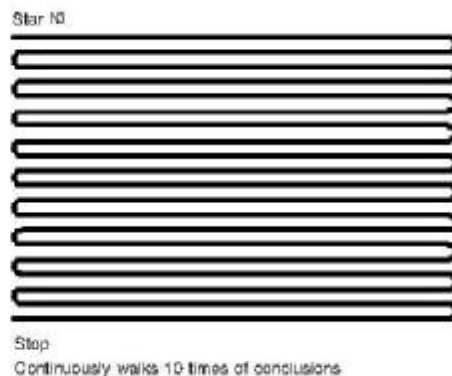

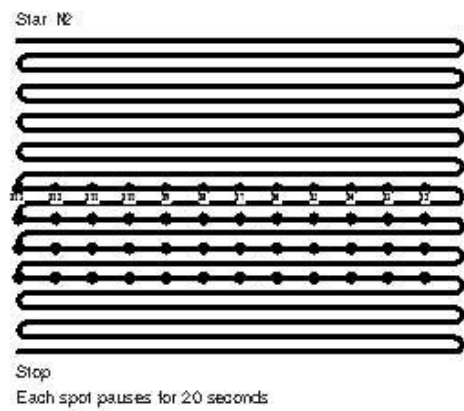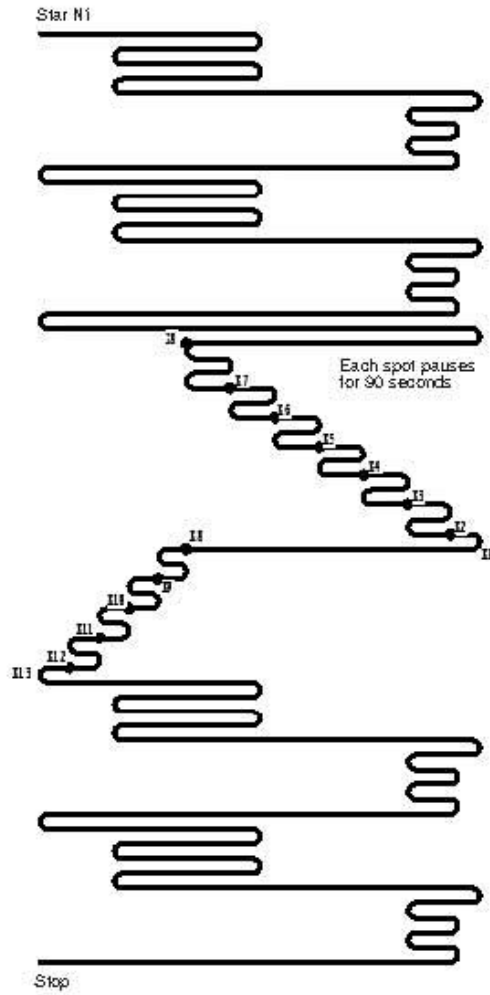


Tension Adjustment




1. Please locate the zipper at the top end of the leather cover.



2. Open the zipper.

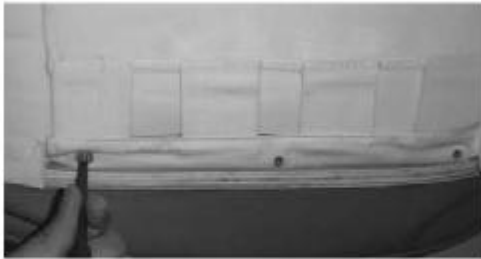

3. Remove all three screws by applying the proper screw driver.

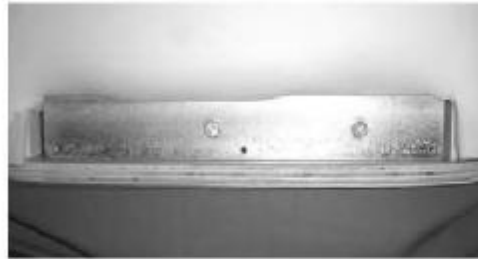

4. Remove the cloth filing as below.



5. Use number 10 wrench to adjust the two nuts as seen in the above picture.



6. Use the wrench to turn the left nut 360 degrees in clockwise.



7. Use the wrench to turn right nut 360 degree in clockwise. the nuts will cause the cable to break over too much tension.

Please not turn the nuts over 540 degree. Over tighten the nuts will cause the cable to break over too much tension.

Installation Of The Product In Your Home

Points To Remember Before Installation

- 1.Product must be placed on the leveled ground
- 2.Please use the product in clean environment
- 3.Please do not place the product near fire
- 4.Please avoid to pour the fluid on the product
- 5.Do not place the product on the bed, couch or any other unstable surfaces

Installation Instruction (See the Set-up procedure)

Maintenance & Cautions

Safety Practice

Safety

- 1.This product has been designed for your maximum comfort and safety.However,improper use may result in fire and electricity shock.
- 2.Please make sure your power output is AC120V~220V before each usage.
- 3.Please use Proleader power cable that shipped with the product. Do not use other power cables other than supplied.
- 4.Please go to your local store or your service provider to replace the power cable if it is damaged.
- 5.Do not overload your power.
- 6.Please do not unzip or open the leather cover in case of danger when the bed is in idle or operation.

Care And Storage

- 1.Please unplug the power cable after each usage.



2. Do not attempt to open the product if you are not a professional service engineer.
3. In the case of abnormal noises, please stop using the production immediately and contact the service people.
4. Please do not mix use all cables to avoid danger.
5. Please do not pull directly on power cable when it is difficult to unplug to avoid danger.
6. Do not use wet hands to touch cables to avoid electricity shock.
7. Use separate outlet with other electronics.
8. When idle for a long period of time or in the storm weather, please unplug the power cables from your power outlet.

Remote Control Operation & Error Messages Explanations

The concrete question please examine the maintenance instructions handbook

Troubleshooting & Service

Self Troubleshooting

Problems	Check	Solution
Remote control has no or wrong display	Power cable plug in properly	Reconnect power cable to outlet
	Remote control plug in properly	Reconnect remote control cable to the outlet
	Burnt fuse?	Replace the fuse
	Power switch on?	Turn power switch to On position
	Press wrong Function buttons	Read the instruction
Projector does not light up	Check if power is plugged in and outlet has power	Call service if none of these work. Do not open the product yourself
Internal projector does not move	Press the remote control to choose the programs	Shut off power and restart the product
	No indicator light	Shut off power and restart the product

Service

Please contact your local store or dealer for service information if problems persist. Please refer to them for the service charges.

Specifications

Product Name	Jade Thermal Massage Bed
Model	
Voltage/Frequency (For Europe)	220V+/- 20V / 50HZ+/-1HZ
Voltage/Frequency (For NA)	120V +/- 20V / 50HZ+/-1HZ
Watts	<= 400W
Working Temperature	Temperature from room temperature to 70°C
Lot No.	See machine plate
Manufacture Date	See machine plate
Specs	See machine plate
Net Weight	
Gross Weight	135kg
Manufacture	
Humidity	

Additional information

Working diagram

