

Thermal Therapy Apparatus

Fan-infrared health products

MANUAL





Features

- 1、 the display is clear and comprehensible.
- 2、 easy-to-use operation panel.
- 3、 temperature settings for units 1 °C, divided into 3, 6, 12-hour time set.
24-hour automatic shutdown.
- 4、 the temperature in 20 °C --- 70 °C range can be set free. With features to prevent overheating, the failure showed that when the temperature [88], and stop heating.
- 5、 with negative ions and electromagnetic Qiaduan function.
- 6、 a unique self-detection.

Technical specifications

- 1、 suitable for 0 °C ~ 50 °C environment use
- 2、 the use of power supply: AC220V ± 10%
- 3、 the use of power frequency: 50 Hz ± 10%
- 4、 the use of power consumption power: 60 W ± 10%

Xiuyu


Xiuyu a serpentine quality jade, Variety of colors, mainly to Danlv, both yellow, white, brown and other colors, More delicate texture, A greasy or shiny wax-like color and translucent. The magnetic can produce excellent biological role, to improve circulation, stimulate regeneration, Increased activity, physical function, enhance phagocytosis and the formation of antibodies. Xiuyan Yu is the human body has an excellent health care role.

Zheyu

Ge-containing the body of selenium, zinc, nickel, cobalt, manganese, magnesium, calcium more than 30 kinds of trace elements helpful to the human body, these elements distribute the start-up wave and the start-up wave of human cells with a volatile state, the human body As the cells from the germanium-distributed to resonate fluctuations and resonance, so that human cells more dynamic organization, and to promote through the bad blood, enhancing metabolism, and body waste into the exclusion. Ge on human medical fitness, "Compendium of Materia Medica - the Department of Epigraphy Vol. VIII" records, this-is: "Apart from stomach heat, Chuanji trouble Men, flexible tendons strong bones, Zhike, Run-heart, and assist the sound jets "Effect. It can produce excellent biological role, to improve circulation, stimulate regeneration, increased activity, physical function, enhance phagocytosis and the formation of antibodies. Therefore it is an excellent human health effects.

Black Lvyu

Lvyu also known as your black serpentine, it rich in a lot of activity and the body of magnet for dozens of trace elements minerals, and the effectiveness of a two-way adjustment, the plot may be in addition to heat, heart and lung Run, spleen and stomach-benefit Soundjets, mayonnaise hair, raising five internal organs, Lee blood, the eyes and ears, long used light-weight for many years, known as the "medicine Wang." This jade in a specific temperature, can release large amounts of the far infrared and physical requirements of various trace elements, can change the quality of body fluids, the body required to achieve the health of PH value of the environment, improve the human immune function integrated, adjustment Eight humansystem, improve sleep quality, on hereditary diseases, chronic consumption of diseases, cardio-cerebral vascular diseases, nervous system and endocrine system and diseases such as with unexpected results. It not only has high ornamental value and collection value, more other-class unparalleled therapeutic health value, so it is extremely valuable and precious stones by the World Association of jade found in the world Need.

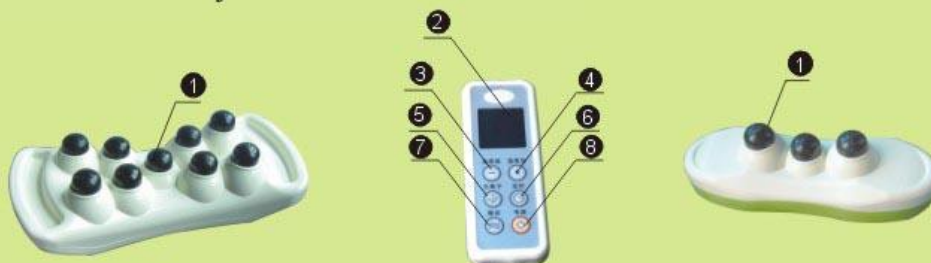

- | | | |
|--------------------------|-------------------|----------------------|
| 1、 3、 9 massage ball | 2、 LCD display | 3、 temperature(down) |
| 4、 temperature(up) | 5、 anion function | 6、 time setting |
| 7、 vibration mode select | 8、 power switch | |

The basic methods of operation

- 1、 Turn on the power: the status is working, setting temperature is 45℃ and the anion function is OFF, vibration is OFF, auto power off function is automatic off after 24h.
- 2、 Adjust Temperature: the temperature is depend on the personal situation and different massage place, you can press the up and down button to adjust it, one press is one degree, and you can see the number on the control panel. Do not stay one place long time.
- 3、 Adjust anion function: press the button, it is on and can see the sign on the panel, press again, it is off and the sign will void.
- 4、 Adjust auto-off function: press each button, the auto-off time will add 10 minutes, long press is working, the max is 720 minutes.
- 5、 Vibration function:
 - a. strong vibration: powerful range, do not use at neck.
 - b. week vibration: lower range, it can use all of the body.
 - c. periodic vibration: 4 second per one cycle, it is switch at ON/OFF, the vibration is strong.
 - d. Stop state
- 6、 Power: press it, the massage is working or standing.
- 7、 Please do not plug the power when not using

Technical specifications

- 1、 suitable for 0 °C ~ 50 °C environment use
- 2、 the use of power supply: AC220V ± 10%
- 3、 the use of power frequency: 50 Hz ± 10%
- 4、 9 ball use of power consumption power: 60 W ± 10%
- 5、 3 ball use of power consumption power: 30 W ± 10%

Far-infrared health products **Warranty cards**

Model:

| User name | Address | | |
|------------------|---------------------|-------------------|-------------------|
| Purchase stores | Invoice number | Date of purchase | |
| Maintenance date | Maintenance records | Replacement parts | Maintenance staff |
| | | | |
| | | | |
| | | | |