

Thermal Therapy Mat

TOURMALINE HOT PAD

PRODUCT MANUAL



Thermal Therapy Mat

TOURMALINE HOT PAD

THE STRUCTURE AND USAGE OF TEMPERATURE REGULATOR



Usage:

1. Turning ⑤ with clockwise direction
① is on, it means the power is put through, turn with anticlockwise direction, it means turning off the power.
2. turn ⑤, select the moderate temperature; when ④ is flashed to set the temperature; after setting, it displays the actual temperature of hot pad.
3. Static electric testing touching key ⑥ if ③ is on, it means it produces static power, touch ③ again and again, it can get rid of electromagnetic wave.

CAUTIONS

- 1) Before operation, read the operation manual carefully, and operate it according to the indication of operation manual.
- 2) Alarming- before operation, check the power, the voltage is qualified with the sign or not.
- 3) When operating, it must ensure the electric heating pad and soft wire without any wrinkle or blocking phenomenon.
- 4) Do not insert pin etc. hard objects.
- 5) To avoid damaging burning person, prohibit using by disabled person, baby or unsuitable person.
- 6) If this product is put on the folded quilt etc. thick and heavy object. Do not put through power to avoid over-heating on some part.
- 7) It is prohibited to use with other heating equipment (such as foot oven, heating bad) together, prohibited to heat for other objects.
- 8) It is prohibited to air, toast and iron badly to avoid damaging the insulation of heating wire.
- 9) No sparkling (such as smoking etc.) to avoid causing fire.
- 10) If it is not used for a long time, cut off the power and take the plug off to avoid accident.
- 11) Do not squeeze the wire or operate the switch by hand, do not pull the tool by it.
- 12) It is not suitable for children to use. Unless parents or supervisor do it well for them. Or guide the children to know how to use it.
- 13) When the electric heater is put through power, do not cover the control device.
- 14) When operating the heater all night, please set it with lower gear.
- 15) This product can only be used with the unloading device harmonized with it.
- 16) Check the heater, especially the power wire, whether it is damaged or not. If the power wire is not good, replace the soft wire produced by our company to make it matched well with each other.
- 17) Prohibit washing up or dry washing, when damping, do not use it.
- 18) When it is not used for a long time, pack it well according to original packaging method, keep it in the cool and dry place, do not squeeze it fold it to avoid heating and bending seriously.
- 19) If something trouble happened (such as smoking, fire, irregular sound or smell etc) cut off the power first, then deal with others.
- 20) In order to make it work normal, and prevent other accidents happening, please obey the regulations above mentioned, otherwise all the accidents must be in charge of yourself.

PRODUCT GUARANTEE CERTIFICATE

Product Guarantee Certificate			
Users' information	Name		Contact telephone
	Add		
	ID numbers		
Product information	Name		Model
	Date of purchasing		Agent




Tourmaline ore



Extract high-quality parts to make powder granular.



Extract pure and light of Tourmalin powder

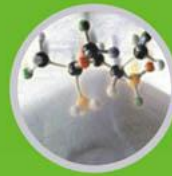

Release a large number of negative ions



Select high-quality of Tourmalin powder to meticulously make production molding.



Under the 1350 high temperature to burn for 48 hours



Make out the glossy Tourmalin after a special heat treatment and release far-infrared and negative ions.



1

Promote metabolism



2

Regulate vegetative nerve



3

Promote blood circulation



4

Regulate endocrine function



5

Reduce blood viscosity



6

Hypolipidemic



7

Improve the function of internal organs



8

Lift fatigue



9

Improve blood circulation



10

Enhance immunity mechanism



11

Neuralgia



12

Rheumatoid arthritis



13

Cervical spondylosis



14

Scapulothoracic periarthritis



15

Back pain



16

Acute and chronic injury of soft tissue



17

Vegetative nerve system dysfunction



18

Gastric wrist pain



19

Prostatitis



20

Pelvic inflammatory disease